

User Manual



Thank you for purchasing the Sonic Cleanser - Blade Style. To use the Sonic Cleanser correctly, please read this user manual before use.



Product Structures & Parameters




Name	Sonic Cleanser - Blade Style
Technology	Ultrasonic at 24,000 pulsations/sec
Working Voltage	3.7V
Output	5V/1A
Rated Power	3W
Mode 1	Cleansing Mode
Mode 2	Treatment Mode

Product Instructions


Cleansing Mode:

1. Apply cleanser on damp skin;
2. Turn on the sonic cleanser, and slowly move the curved blade probe downwards along the facial skin. The dirt and grime on skin are purged by high-frequency ultrasonic wave.
3. After use, turn off the device and clean the sonic cleanser blade probe with a clean cloth.

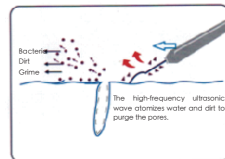

Treatment Mode:

1. Liberally apply Mineral Booster to the face and neck;
2. Apply sonic cleanser curved blade probe upwards for antioxidant infusion and turn on the device. The inverted probe contacts the skin at an angle of 45 degrees. Slowly glide the blade probe along facial skin to enhance skin absorption and stimulate skin radiance;
3. After use, turn off the sonic cleanser and clean the device with a clean cloth.

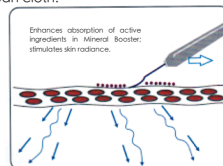

Product Usage Notes:

1. The sonic cleanser blade probe should be used at an angle of about 45° to the skin surface;
2. Avoid sensitive areas such as the throat, eye and mucosa, etc.;
3. The blade probe should be used with continuous movement along the skin surface and should not be placed stationary at a single point when being used.
4. Clean the blade probe with a moist cotton pad before use.
5. The device should be fully charged before first use. Indicator is lit red when charging and green when fully charged; the device will take approximately 3 hours to fully charge.
6. The sonic cleanser has a fixed usage time. The device will automatically turn off after about 10 minutes of use. Please restart the device if needed.

Safety Precautions:

1. Please turn off the device before cleaning the product;
2. Only the sonic cleanser blade probe can be washed with water, do not immerse entire device in water;
3. Do not expose the device to high temperatures or heat as it has a rechargeable battery;
4. Please do not use this device if you are allergic to metallic products;
5. Please stop using the device if it becomes warm or hot;

Do not use this product in the following cases:

- A. During pregnancy or postpartum lactation;
- B. Allergies, eczema, skin diseases, skin inflammation and sunburn;
- C. Any form of open wounds or skin damage.;
- D. Post-surgery or recuperation after illness;
- E. When feeling unwell or running a fever;
- F. When drinking alcohol or taking medications.